
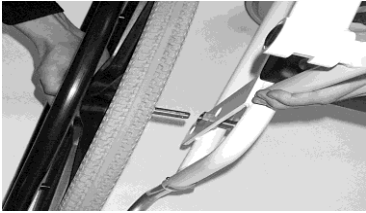
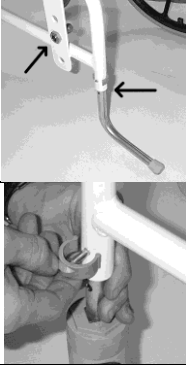


User guide	
M2 shower/commode chair with 24 inch wheels	
How to mount the wheels:	
	The axel is clicked into the bush. Pull the wheel to make sure it is correctly mounted. But be careful not to tighten the bolts too much. Using the chair without the driving belt, the wheels will get more elastic. (when you receive the chair it is adjusted to the middle position in height)
How to adjust the height:	
The bushings to the 24-inch wheels are taken off with a wrench key and moved to the desired height. The anti-tip bar is moved to the same position as the 24-inch wheels. The clips holding the front wheels are taken out, and the front wheels are moved to the desired height. The front wheels can be pu in a position which is ONE (1) position different from the 25 inch wheels (if you want the chair to be slightly tilted) All other adjustments shall be in the same position! (For example: in the top hole). Double check that all nuts and bolts are tightened and that all the clips are pressed in place.	
Make sure to tighten the bolts well	
	With the help of an unbraco key you tighten the screws so they give the wanted brake impact. Tighten the bolts well!
Footrests	
Pushing the footrests aside: Lift them and turn them to the side you want. The plastic plate can be tipped up if you want to take them to take up very little space. Height adjusting the footrests: Adjust the height of the footrests by taking the clips out, placing the footrests in the wished position, and last put the clips back again.	
Pan / bucket	
	When sliding the bucket of the chair, it is important to slide it straight backwards, along the seat. Otherwise it can be very difficult to get the bucket of the chair.

WARNING

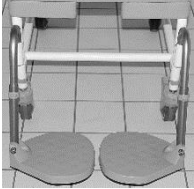

RISK OF ROLLING OR GLIDE-OUT!

TO AVOID THE CHAIR CAN ROLL OR GLIDE IN TRANSFER,
TO ALL 4 WHEELS BE LOCKED!

ATTENTION !!

IF THERE IS SMOOTH SURFACE,
THERE IS STILL A RISK OF THE CHAIR CAN ROLL OR GLIDE.



DO NOT SIT ON THE OUTER 10
CM OF THE SEAT THAT THE
CHAIR CAN TIP FORWARD.



DO NOT STAND ON THE
FOOTRESTS



Check with proper intervals the fastening of the wheels to the leg.

2 Years complaint rights against manufacturing errors.

Terms of claim for frame.

HMN grants 5 years claim rights against frame break-down for all M2
Shower/commode chairs and stainless steel shower stools.

The warranty applies for production and material break downs, and can
only be used if the products have been in normal use.

Please note that the normal wear is not included in the warranty.

Cleaning:

The product can be cleaned with regular cleaning without
non-abrasive detergent. Wash and dry the product after the
cleaning. The product can be washed up to 85° degrees in
3 minutes.

11.05.15 – produced in the EU



HMN a/s

Midtager 22, 2605 Brøndby, Denmark

Tel.: +45 43 96 24 00 Fax. +45 43 96 26 92 E-mail info@hmn.dk

www.hmn.dk